





# INTERDENTAL- BÜRSTEN COMFORT

GB Comfort interdental  
brushes

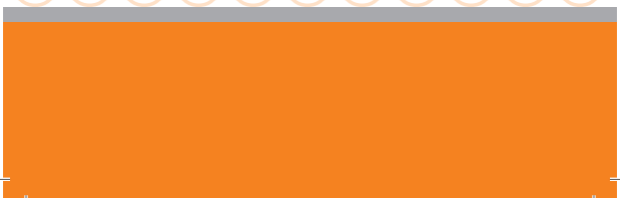

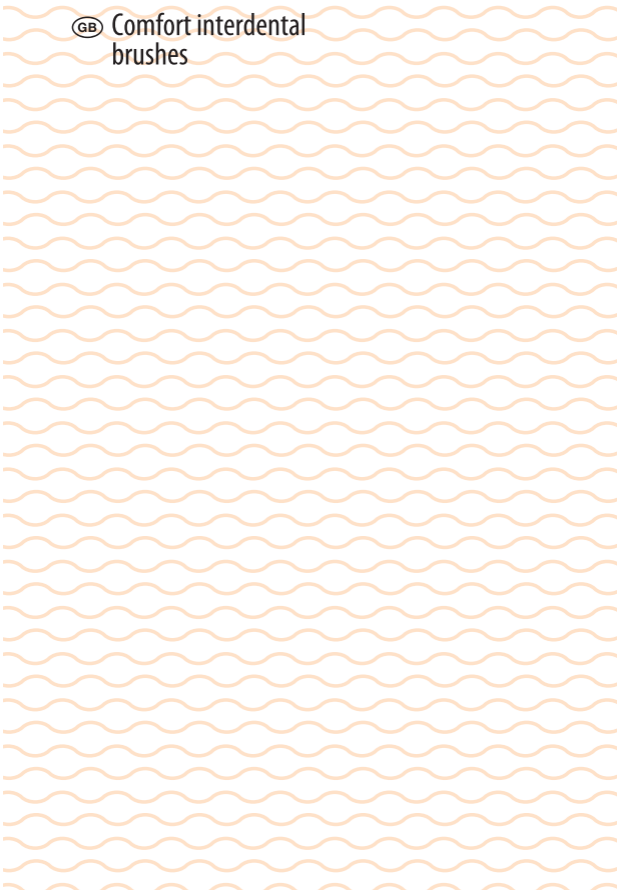

# INTERDENTAL- BÜRSTEN COMFORT

## Anwendung

1. Führen Sie die Interdentalbürste vorsichtig zwischen die Zahnzwischenräume.
2. Bei der Wahl der Interdentalbürste ist darauf zu achten, dass die weichen Borsten die Zahnzwischenräume ausfüllen ohne den Zahn mit dem Drahtkern zu berühren.
3. Bewegen Sie die Interdentalbürsten vorsichtig in den Zahnzwischenräumen zwei-bis dreimal hin und her um Plaque und Essensreste zu entfernen.
4. Nach dem Gebrauch empfehlen wir die Interdentalbürsten mit Wasser abzuspülen und zu reinigen.
5. Benutzen Sie eine neue Interdentalbürste sobald sie Spuren von Abnutzung zeigt.

## GB COMFORT INTERDENTAL BRUSHES

### Direction for use

1. Insert the Interdental brush with care in the wider spaces between teeth.
2. When choosing the size of the interdental brush, make sure that the soft bristles fill out the interdental space without touching the teeth with the wire.
3. Use a gentle, straight, back and forth motion 2-3 times to remove any plaque or food remnants.
4. After use we recommend rinsing with water.
5. Use a new brush when the bristles show signs of wear.

Size	Colour	Diameter Brush	Diameter Wire	REF
superfine	white	1,6 mm	0,40 mm	000 889
xxxxfine	pink	1,8 mm	0,42 mm	000 881
xxxfine	orange	2,0 mm	0,44 mm	000 882
xxfine	red	2,5 mm	0,45 mm	000 883
xfine	blue	3,0mm	0,45 mm	000 884
fine	yellow	3,5 mm	0,57 mm	000 885
medium	green	5,0 mm	0,57 mm	000 886
large	purple	6,0 mm	0,77 mm	000 887
xlarge	grey	9,0 mm	1,06 mm	000 888

Stand: 140717

**M+W Dental**

Müller & Weygandt GmbH • Reichardsweide 40  
D-63654 Büdingen • Tel.: +49 (0)6042 – 88 00 88